



OVERLINE

ITALIA



Infinity

Made In Italy  Technology 

LYMPHO SYNERGY

STARTING OVER WITH ENTHUSIASM

-  PRESSOMASSAGE
-  LYMPHATIC DRAINAGE
-  ELECTRO REMODELING
-  WARM SYSTEM



Dimensions: L600 x P570 x H1260

Tone, lightness and vitality

Overline Infinity Lympo Synergy is the innovative Overline equipment for slimming the silhouette, which uses the synergy of heat and Presso-Massage enhancing the benefits and aesthetic effects of both. A unique system, a great novelty for the aesthetic sector. The effects of Overline Presso-Massage are widely known. Its leg and arm progressing inflation cuffs and layered inner air chambers prevent the veno-lymphatic return, simulating the passage of the hands and ensuring a uniform massage. The Overline innovation lies in combining Presso-Massage action and heat reducing action, a happy union that enhances the treatment, making it more effective and more pleasant.

The final result of Overline Infinity Lympo Synergy is the progressive decrease in circumferences, thanks to the resorption of the oedematous state and to the improvement of veno-lymphatic circulation. The iSol muds and bandages, designed ad hoc for treatments with Overline equipment, play an important complementary action.

Included

3130-12-PWF-INF 12 SECTORS DE LUX model leg cuff

(Legs, Abdomen & Arms)

3130-08-PWF-INF 8 SECTORS DE LUX model leg cuff


(Legs & Abdomen)





SUPER SILENT TECHNOLOGY

This equipment is endowed with a new generation engine that makes the treatment quieter and your work efficient over time with a duration of 5000 hours equal to 10000 treatments.

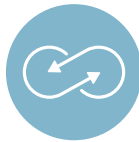

6 PROG. BODY - GENERAL

Treatments that allow you to work simultaneously on the lower and upper limbs. Diversified programs, according to individual needs and customer issues.

- 1) PERISTALTIC CLASSIC SOFT
- 2) PERISTALTIC CLASSIC INTENSE
- 3) MIXED SOFT
- 4) MIXED INTENSE
- 5) MIXED UPPER
- 6) MIXED LOWER

4 PROG. ARMS


Ascending peristaltic treatment 4 sectors.
Promotes the discharge of liquids, reshaping the part treated.

- 1) PERISTALTIC CLASSIC SOFT
- 2) PERISTALTIC CLASSIC INTENSE
- 3) SWOLLEN / HEAVY ARMS
- 4) SPECIAL DREN

12 PROG. LEGS


8 sectors treatments that act on the lower limbs.
Diversified programmes, according to individual needs and customer issues.

- 1) PERISTALTIC CLASSIC SOFT
- 2) PERISTALTIC CLASSIC INTENSE
- 3) CULOTTE DE CHEVAL COMPACT
- 4) CULOTTE DEC HEVAL SOFT
- 5) COMPACT DREN CELL
- 6) SOFT DREN CELL
- 7) EDEMA DREN (40MMHG) VODDER
- 8) ABDOMEN
- 9) PREGNANCY (on prescription)
- 10) INTENSE DREN LYMPH
- 11) SOFT DREN LYMPH
- 12) EDEMA DREN LEGS (40mmHg)

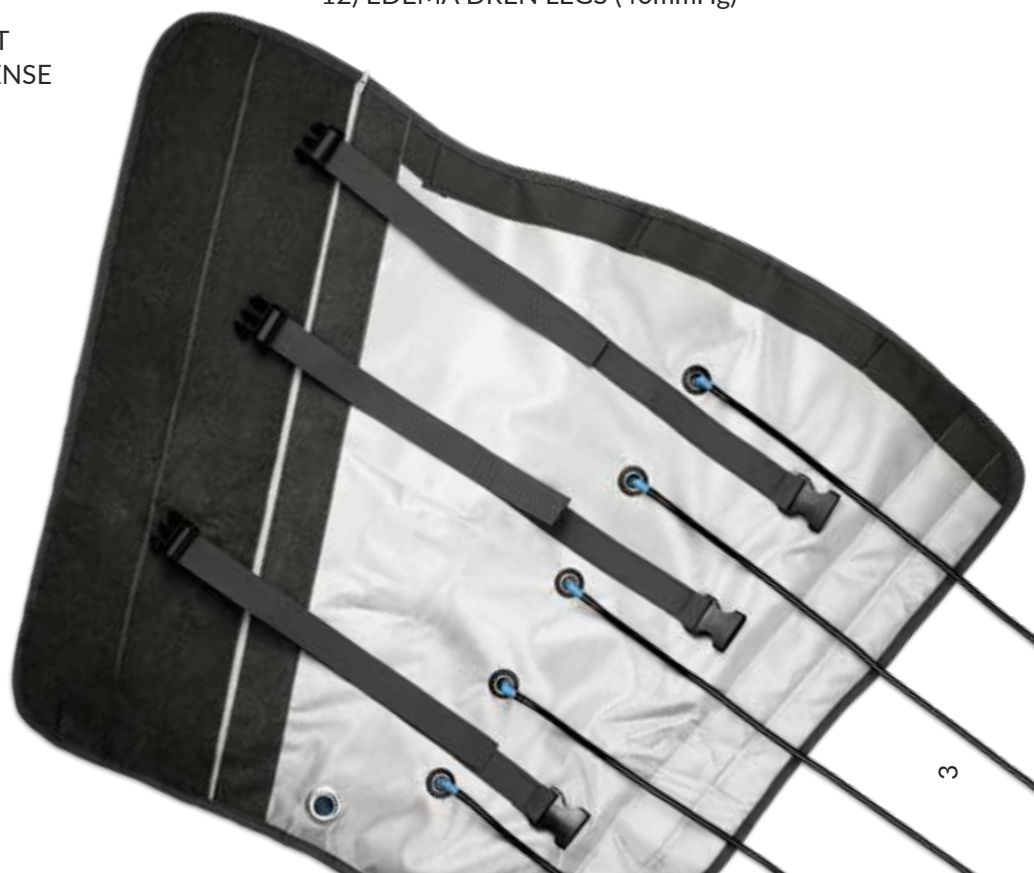

PRESS MASSAGE

Revitalizing Action

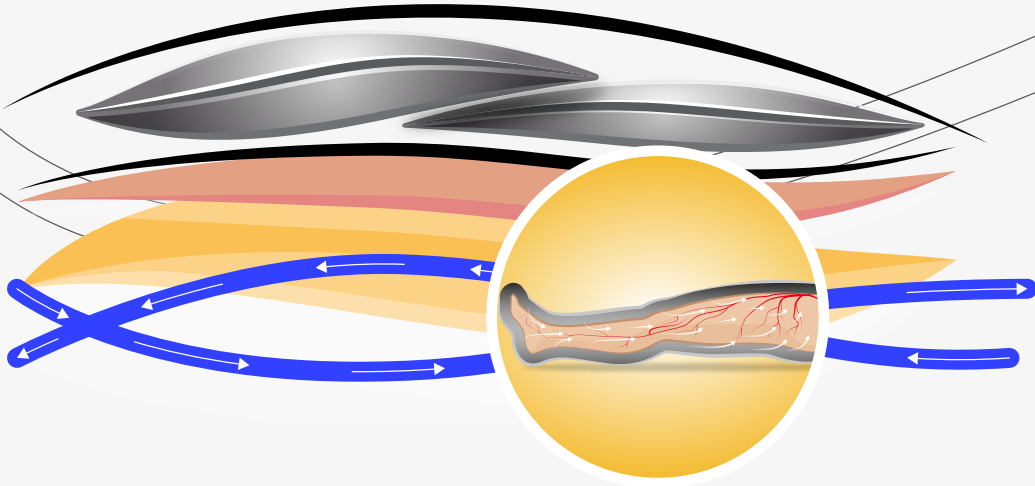


The Presso-Massage contributes to the restoration of veno-lymphatic function, promoting peripheral circulation. It has a draining activity, reducing swelling and muscle tension. It has an important effect detoxifying, revitalizing, improves the oxygenation of tissues and prevents the formation of imperfections caused by cellulite. This treatment is particularly relaxing and generates a sense of well-being and lightness.



PULSE MASSAGE

Draining Action

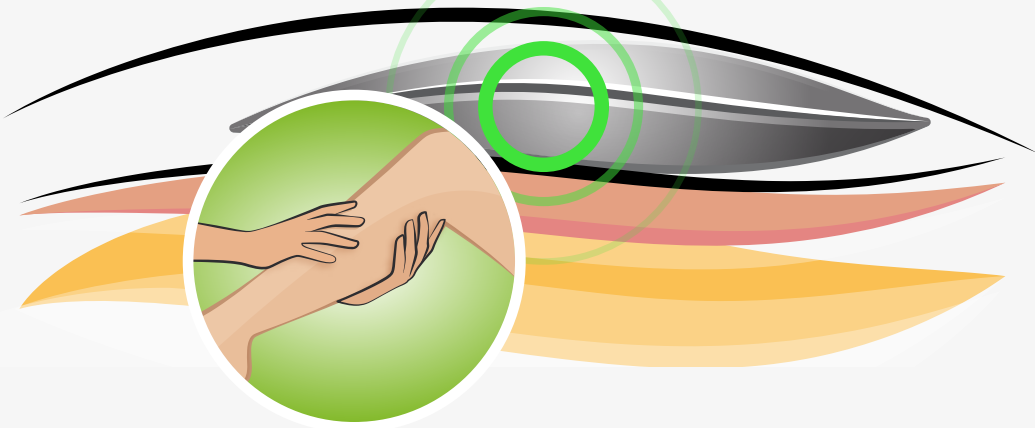


The new PULSE system of Overline Infinity Lympho Synergy was born from the studies of the Overline scientific laboratories in order to make the drainage/massage treatment carried out by the equipment even more pleasant and more similar to manual massage.

PULSE technology uses a computerized system managed by Overline Infinity Lympho Synergy able to perform a pulsed and continuous pressure in every sector. The pulsed and intermittent pressure creates a greater stress, improving the effects of Presso-Massage treatment. The stimulating action on the subcutaneous tissue improves contractility and muscle metabolism, acting on the reflex route on the arterial circulation and activating a natural drainage mechanism.

PRESS WARM


Slimming Action

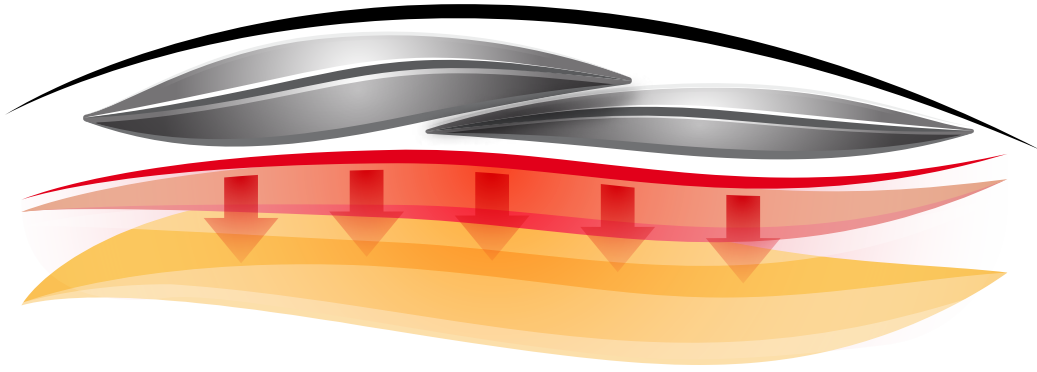

The heat generates pleasant effects and real benefits on the body such as: increased microcirculation resulting in improved interstitial exchanges and increased local metabolism. The increased metabolism has a trophic action, that is, a tissue nourishing action, but also of oxygenation, an activity that acts effectively on the imperfections caused by cellulite. The mechanism of exchange and activation of (local) metabolism is able to promote detoxification. Finally, the heat has important activation properties of the active ingredients contained in the Active Muds and the specific Bandages of iSol Beauty.

DEEP PHYSIOSTIMULATION

Toning Action

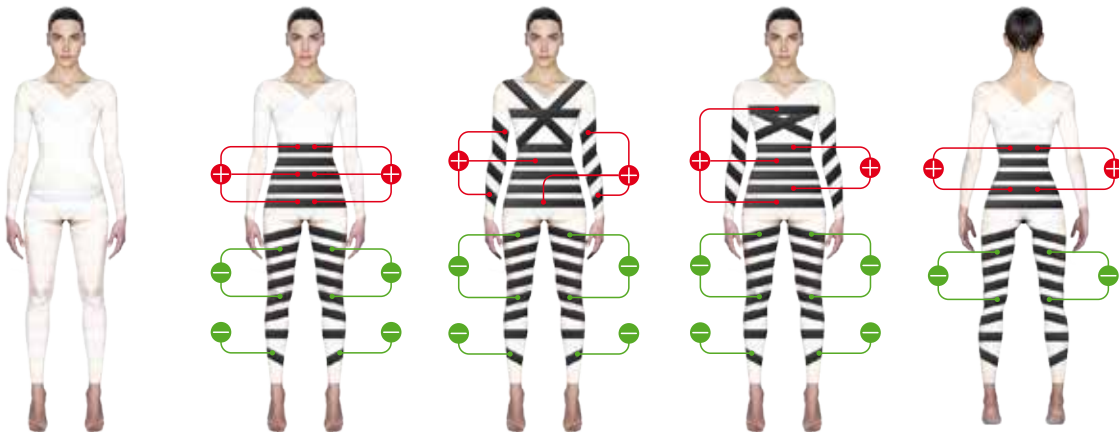

The Modeling Physiostimulation uses frequencies ranging from 30 to 80 Hz. Ideal for toning and recommended in the presence of sagging tissues, this treatment acts on the entire length of the muscles, keeps the skin supple, stimulates the flow of lymph improving its transit, further speeds up the flow helping to dissolve the occurring congestionst.

Remodeling treatment: preparatory steps

The bands for Physiostimulation are to be applied above the iSol Beauty bandages which have moisturizing, toning, anti-cellulite, reducing, draining and anti-fatigue properties.

EASYLOCK

CONNECTION SYSTEM



One of the most common problems in the connection of cables for the at-massage is their insertion not perfect. In this case the technology can not work, causing inconvenience to the customer ready for treatment. Overline, who has always listened to the opinions and opinions of Beauticians, has created "EASYLOCK", a simple but cutting-edge system, which groups all the cables in one connector with special features, which must simply be connected in the socket on the back of the equipment. The tangle of wires will be just a bad memory.

Overline is on the side of the professional Beautician and listens to facilitate her future!



Air chamber positioned under the Lejars's sole and vault triangle to improve the return force.



Leg cuff for lower limbs adaptable both in height and in circumference, thanks to the velcro that offers a perfect seal.




Adjustable triangle to improve pressure at the level of the inguinal ganglia and to manage the proportion to the height of the hips.

It's divided into foot, leg, panty with lymphatic ganglia opening, side-adjustable bodice and includes the convenient pouch to store all the parts that make up the cuffs. The arms sectors are divided into 4 and have a specific chamber for the lymphatic ganglia under the armpit. Adjustable in height and perimeter, and with interchangeable air chamber



Safety belts to ensure perfect adherence along all the areas treated during the expansion of the air chambers

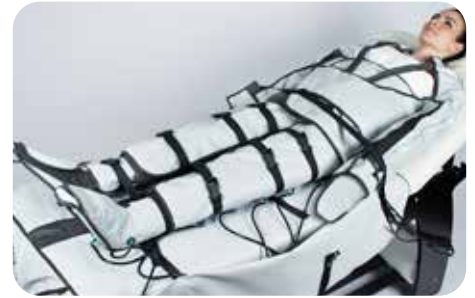

Protective sheet designed with biocompatible materials, washable in complete safety. designed with biocompatible materials, washable in complete safety.



Layered air chambers can be inspected and replaced individually.





VIDEO EDUCATIONAL

The technological heart of the equipment, thanks to the advanced software, offers the possibility to view some videos that help to clarify doubts and improve manual skills. At the forefront of service, Overline Infinity Lympo Synergy is always at the side of the professional aesthetic operator to offer a result of excellence. The video tutorials related to individual protocols and specific manual skills, are easily accessible via the 10" touchscreen. Overline has thought of this solution to allow anyone, in case of a small forgetfulness, or when new operators are involved, to carry out self-training, in order to achieve the skills needed to offer and ensure the customer the best service.



INTUITIVE SYSTEM

Protocols and performances carried out within the Company* respecting the suggested working methods, confirm the synergy between Overline technology and iSol.Beauty cosmetics, ensuring maximum treatment effectiveness.




SISTEMA INTUITIVO

The advanced software of Overline Infinity Lympho Synergy allows you to enjoy numerous and advanced functions, and to operate the equipment via the 10" monitor in a simple, fast and intuitive way. This monitor, as conceived, immediately becomes a fundamental tool to implement the work of the professional.

With a few small gestures, the operator is able to carry out a 360° consultation of the specific and physical situation of the customer, can examine the images, can control programs and treatments, can have access to customer data, and also can verify the situation of the treatment steps to be performed, having the opportunity to motivate and involve the customer more. It is also possible, in the guided path, to watch some videos that help clarify doubts and improve dexterity.



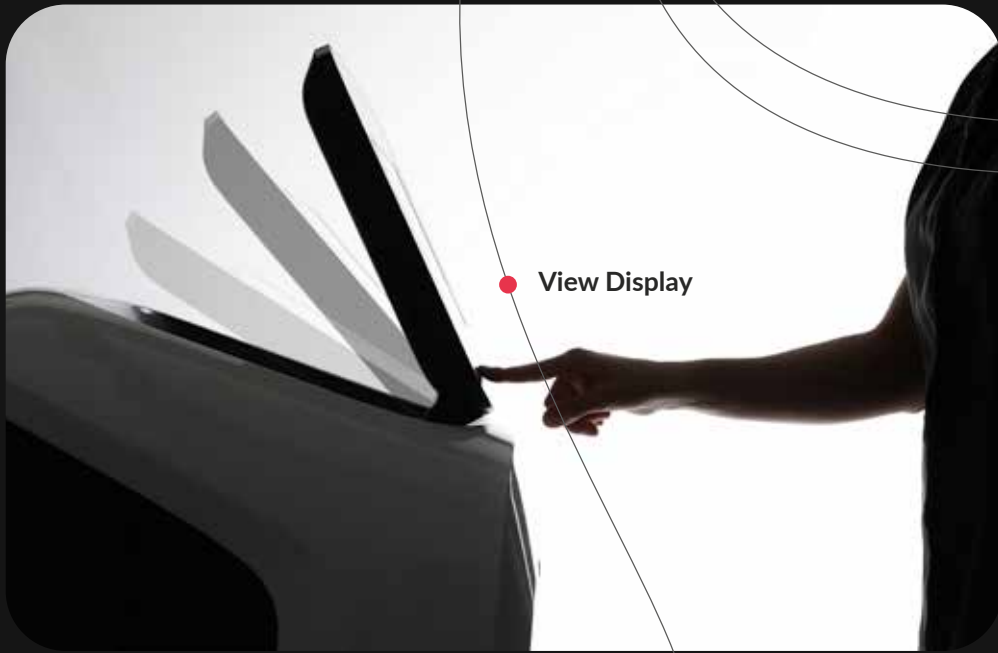

INFINITELY *Beautiful*

Overline Infinity is the innovation choice for the beauty Salon. It is the possibility of infinite customised treatments with the synergy of latest generation technologies. It features an intuitive interface and simple operation. It means effective treatments and Salon customer satisfaction. It involves on-going training support allowing every beautician to make full use of all of the potential offered by the equipment. It provides the safety of using strictly tested, certified systems built in compliance with regulations in force. It has an appealing design of sinuous lines that suits every salon.

Infinity is for those who never lose sight of their end goals!

THE DETAILS

The electronically reclining monitor of Overline Infinity allows easy programming and control of the work, offering a clear vision of the treatment to the customer. A large upper drawer and a storage compartment allow you to store materials and cosmetics, while the shaped shelf on the back of the equipment allows the tidy arrangement of the handpieces.



Storage compartment



Large top drawer

The customization concept for Infinity Overline is very important. Each technology has the advanced Contactless Card reading system, with which, for each customer, the treatment sequence is inserted, storing the steps already made and those to be performed, optimizing the time.



Contactless Card





GARANZIA 36 MESI
ZERO
PENSIERI



The Overline *Infinity Limpho Synergy* is guaranteed from material and workmanship defects, ex-factory, for a period of 36 months. Activate the conditions of validity of the guarantee is very simple, just fill in all its parts the Warranty form that is inside the package of *Infinity Limpho Synergy*, and send it to Overline no later than 10 days from the date of purchase.



INTEGRATED COSMETICS

The combination Overline Infinity Lympho Synergy and the integrated cosmetics of isol.Beauty is synonymous with concrete answers to all aesthetic needs. isol.Beauty cosmetics are designed specifically to work efficiently with the technology of which the equipment is endowed, for maximum operator satisfaction and an excellent beauty experience for the customer.



1 2 3 4 5 6 7 8

- | | | | |
|--|---|--|--|
| <p>Muds</p> <p>1 FANGO-SLIM
Cream mud with white clay, fenugreek and dead sea salts
ISO.MUD.200
[1000ml]</p> | <p>Muds</p> <p>2 FANGO-CELL
Cream mud with green clay, active algae and green tea
ISO.MUD.100
[1000ml]</p> | <p>Conductive gel</p> <p>3 DRENA-TONE GEL
Conductive gel toning and reactivating effect, with ginseng
ISO.TONE.200
[1000ml]</p> | <p>Conductive gel</p> <p>4 LIPO-CELL GEL
Conductive gel for cellulite and fat, with guarana and green tea
ISO.DREN.100
[1000ml]</p> |
| <p>Cleansing</p> <p>5 CLEANSING SOFT MOUSSE
Skin Compatible Face & Body Cleanser
ISO.FOAM.150
[150ml]</p> | <p>Emulsion</p> <p>6 HYDRABODY CELL
Body emulsion for cellulite
ISO.BODY.100
[200ml]</p> | <p>Emulsion</p> <p>7 HYDRABODY FIT
Body emulsion toning effect
ISO.BODY.200
[200ml]</p> | <p>Scrub</p> <p>8 BODY SCRUB
Salt scrub with argan pits and natural oils
ISO.SIL.002
[500ml]</p> |



9 10 11 12 13 14

- | | | | |
|--|---|---|---|
| <p>Bandages</p> <p>9 DRAIN ACTIVE
Draining Detox
ISO.BE.115</p> | <p>Bandages</p> <p>10 REDU CELL
Cellulite reducing effect*
ISO.BE.125</p> | <p>Bandages</p> <p>11 SLIM REDUCING
Slimming Effect
ISO.BE.120</p> | <p>Bandages</p> <p>12 TONE UP
Volume effect
ISO.BE.110</p> |
| <p>Bandages</p> <p>13 COMPACT FRESH
Treatment with toning effect
ISO.BE.105</p> | <p>Bandages</p> <p>14 THERMO CELL
Reducing effect thermal treatment
ISO.BE.100</p> | | |


* Imperfections caused by cellulite.

Always at the service of beauticians

Cutting-edge and Innovation are the ideal terms to define **Overline**. A historical Milanese company, undisputed leader in the production and sale of hi-tech professional cosmetic equipment.

Founded in 1989 and, ever since, it can very simply be described as a “success story”. A company that expresses its business style through ethical commitment and careful and responsible commercial management.

The literature of the company uses terminology allowed at an international level, as well as being the most commonly used terms to define technologies and features

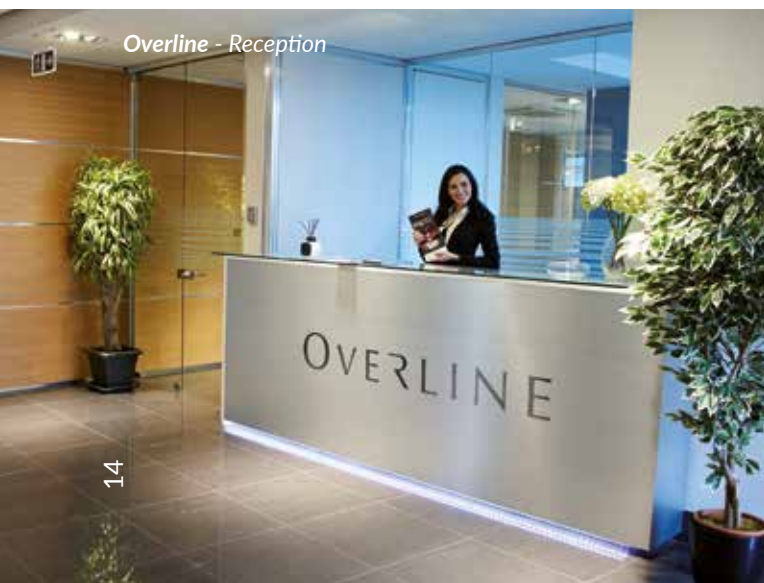

Overline - Reception



Overline - Training course

Research & Development

The beating heart of the company, the innovative Research and Development laboratory, is the place where Technologies of high performance and simple use provide effective solutions to the demands of professional beauticians. Various, safer systems that offer multiple functions for face, body and hair removal

Made in Italy

100% Made in Italy production that offers competitive advantages in terms of performance, safety and service to the customer. Overline technologies are distributed at the finest cosmetic and beauty salons across Italy and in every European country in addition to South America, Mexico, Africa and Asia.

Certifications

Overline is a trademark of Technology S.r.l., a company that constantly pursues the research and development of technologies, to offer the professional beauty market innovative systems that offer visible results, while maintaining the maximum safety of the beautician and customer. Thanks to its great experience, **Technology** possesses important certifications for compliance with regulations and international standards in the production of biomedical, electro-cosmetic and cosmetic equipment. The certifications that make all of the technologies very high level are the following:



Tüv Rheinland ISO 9001:2015. This represents the internationally recognised Quality Management System certification that ensures that the devices comply with the management system required by the standard in terms of design, production, sale, development, and research.

Tüv Rheinland ISO 13485:2016. The certification related to the Quality Management System of Medical Devices is a standard reference for Quality as part of the regulatory procedures governing the design, production and marketing of Medical Devices.

Tüv Rheinland ISO 45001:2018. Relative to the standard for the management of the health and safety of workers.

www.tuv.com ID 9105082302



OVERLINE

ITALIA



Infinity

ARBIDORA

MEDICAL TRADING CO. L.L.C

عربيدرا لتجارة المعدات الطبية ذ.م.م



CALL +971 50 480 7988

www.arbidora.com



Overline reserves the right to make changes to the equipment described, at any time and without prior notice.
Overline equipment is designed and programmed to obtain results on imperfections without exceeding the allowed cosmetic limit.



ISO 9001:2015



ISO 13485:2016



ISO 45001:2018

Overline Infinity is a Technology S.r.l. trademark